

GF. glove factory

GF. cosmetic glove


TEXTILE COSMETIC GLOVES FOR HAND PROSTHESES



GF glove factory GmbH | Breite Straße 155 | 76135 Karlsruhe | Germany

Phone: +49 721 161 784 62 | info@glovesfactory.de

www.glovesfactory.de

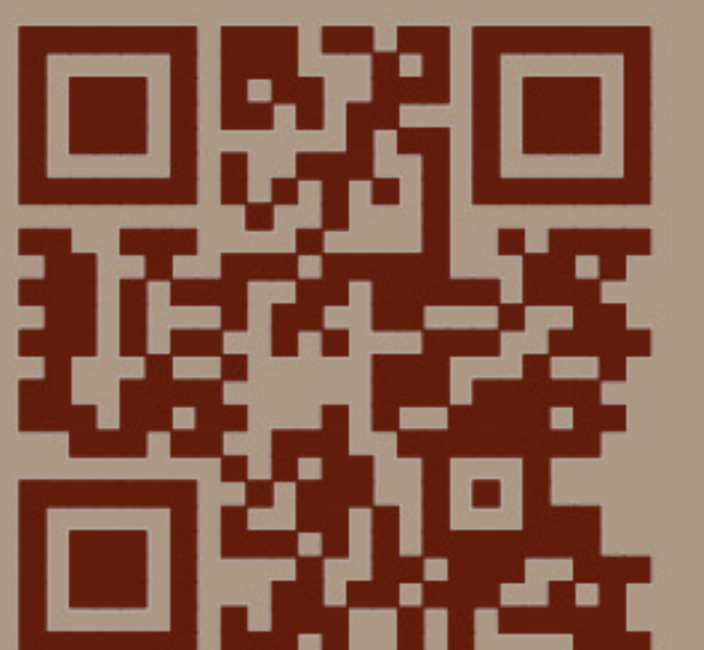

The textile cosmetic gloves **GF. cosmetic gloves** are made of highly elastic knitted fabric with polyamide and Lycra® fibres. Therefore, its especially suitable for high agile hand prostheses with single finger movement. The two-piece glove set consists of a protective underglove and a cosmetic overglove, which reproduces a detailed and natural appearance. This product line of **GF. cosmetic gloves** is specially tailored for the dimensions of VINCENT hand prostheses. Furthermore, special designs for partial hand systems and prostheses of other manufacturers are also possible.

The properties

- Detailed, natural skin appearance with embossed fingernails.
- High elasticity due to polyamide and Lycra® fibres and elastic seams.
- Slip-resistant palm enables safe and non-slip grasping.
- Protection against splash water and dirt due to water-repellent impregnation.
- The overgloves are available in three lengths and different arm circumference.

The advantages

- The high elasticity of the textile cosmetic glove enables unrestricted movements of the prosthesis.
- Minimal wrinkling, as the textile is pre-tensioned to the prosthesis.
- Skin-like appearance.
- Adjustment of the glove colour to your own skin colour by selection of different skin tones.
- Easy glove change.
- Washable.

Additional options

- Individual adaptation of the gloves for partial hand prostheses.

Storage life

- Depending on the degree of use, a cosmetic glove can be worn for 2-8 weeks.
- The durability of a glove package (4x cosmetic glove + 2x underglove) is similar to that of a standard silicone glove.

